


Assembling and cleaning the Pari PEP and Pari LC Sprint

Information for
patients and
parents

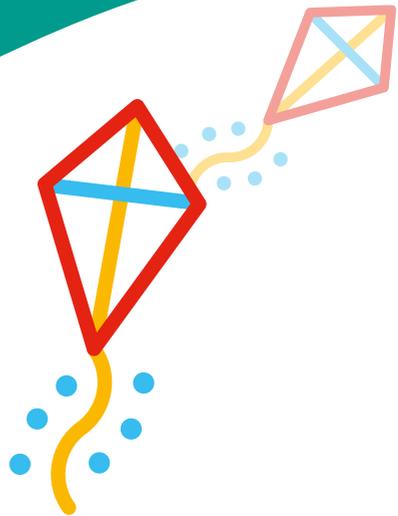

leeds children's
hospital

caring about children


This leaflet explains how to assemble your combined Pari PEP and Pari LC Sprint and how to look after them.



1 Pari PEP mouthpiece




2 Body Section with manometer arm

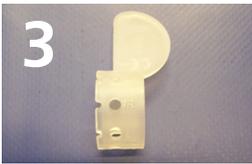

3 Variable hole section



4 LC Sprint top section

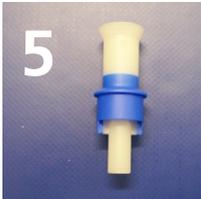

5 Baffle



6 LC Sprint nebuliser pot




7 Compressor



8 Extension tube



9 Pressure gauge

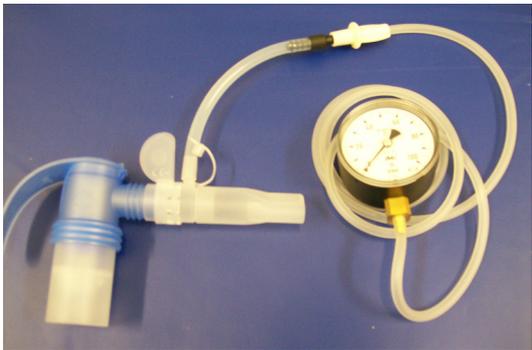

Completed PEP and LC Sprint

Connect the LC sprint to the compressor machine using the nebuliser tubing



Nebuliser tubing
(colour may vary)



Compressor



Please note: in this set up the Pari PEP mouthpiece is used (this mouthpiece has no flap).

Cleaning and maintenance

Your Pari PEP and LC Sprint should be cleaned after every time you use it.

1. Clean parts 1-8 in a clean bowl with warm tap water and washing up liquid. Rinse thoroughly with warm running tap water.
2. Dry all the parts with kitchen roll. Leave to dry in a clean and warm place.
3. Keep the LC sprint parts somewhere clean and dry such as a lunchbox. Put the LC sprint back together just before you use it.
4. Once a week your Pari PEP and LC sprint should be sterilised.

Wash the sections as in Step 1. Then put pieces 1-6 into a pan of water and boil for 10 minutes. **Do not boil parts 7 and 8.**

Carefully remove the pieces from the water and dry as above.

Please bring the Pari PEP and LC sprint to clinic with you each time you attend so we can check if the equipment or settings need to be changed.

Do not wash the nebuliser tubing. To clean this, leave one end connected to the compressor and the other end disconnected, turn on the compressor for two minutes to blow air through the tube. The nebuliser tubing should be replaced every 6-8 weeks. Please ask your physiotherapist/nurse if you need more tubing.

The surface of the compressor and the pressure gauge can be wiped over with a damp, soapy cloth. **Do not submerge the compressor in water.**

The filter in the compressor should be checked every month and changed if discoloured or soiled. It should be changed as a minimum every six months. If you have a different type of compressor your physiotherapist or nurse will give you advice related to the type of compressor machine you have. Please contact your physiotherapist or nurse if you need extra filters.



If you have any questions please contact:

Children's Physiotherapy Team Tel: **0113 392 6361**